

Physiotherapy Devices / Traction Therapy

ORTHOTRAC LASEDIA

Traction Unit / OL-6600

OG Wellness

Feature 1

A simple space-saving design
Easier for patients to get into supine position
Cervical & Lumbar traction unit


Pictures



Pictures

Special Feature 1

Size Comparison

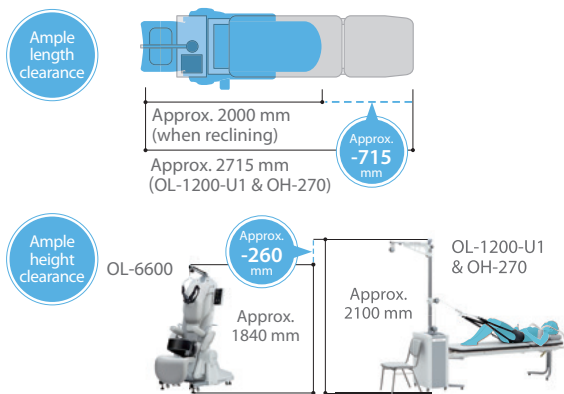

Special Feature 2



Special Feature 3



Special Feature 1

Space saving 2-in-1 design

Perform both cervical and lumbar traction with the patient in a relaxed position, with one "2-in-1" traction unit.

Special Feature 2

A Free Positionable Operation Panel

The control panel can be placed on either the right- or left-hand side upon installation.

Special Feature 3

Coldness Reducing Heater

Equipped with a heater to prevent the patient from feeling a cold sensation around the waist, neck and shoulder areas. Choose between two temperature settings.

Feature 2


Auto-fit side supports
Easier on the patient than seated traction
Cervical traction in a relaxed, supine position

Pictures


Special Feature 1



Special Feature 2

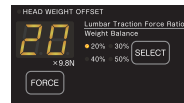

Special Feature 3

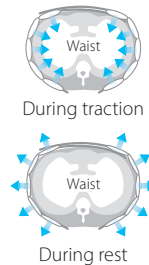


Special Feature 4




Special Feature 5




Special Feature 6

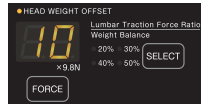

Special Feature 7

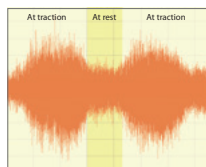

OL-2200 (Seated Traction)

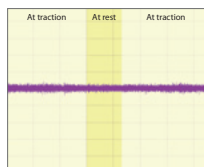

OL-6600 (Supine Traction)

Special Feature 1

Automatically adjusting side supports fitting various body types

The side supports move prior to treating the patient, for a perfect fit. The side supports pull in towards the body to match the traction, providing solid upper body support.

Special Feature 2

Automatic traction force setting

The machine weighs the patient and automatically sets the traction force to the body weight percentage preset by the therapist. This saves the therapist the trouble of calculating the traction force.

- * Can be set manually.
- * Cannot be used as a body scale. The amount of traction force is reduced slightly from the calculated amount for increased safety.

Special Feature 3

Easy-fit inflatable traction belt

Gentle yet firm immobilization of the waist with an inflatable belt. Choose from three tightness settings. In "Relax Mode", the belt automatically deflates during the rest time, to reduce strain on the patient's waist.

Special Feature 4

Easily attachable belt saves therapist's time

The inflatable traction belt uses Velcro for easy attachment. This reduces the amount of time needed to immobilize the patient, so the treatment can start right away.

Special Feature 5

Firm-fit cervical sling

Newly developed cervical sling for a firm fit to the back of the head offering comfortable traction.

Special Feature 6

Head weight offset function

A preset head weight is subtracted from the traction force set for cervical traction.

Special Feature 7

Patients have less muscle tension in the neck in a supine position than in a seated position, resulting in more relaxed cervical traction.

With OL-6600, cervical traction is performed with the patient in a supine position. Compared to seated traction which results in more muscle activity in the patient's neck, a supine position results in less muscle activity in the neck, for more relaxed traction.

Specifications**ORTHOTRAC LASEDIA****OL-6600**

Product name	Traction Unit OL-6600
Rated power supply	AC 110 V/220 V/230 V/240 V 50/60 Hz
Power input	200 VA
External dimensions	Approx. 800 × 1240 × 1840 (H) mm
Weight	Approx. 124 kg
Traction force	Cervical: 29–245 N (3–25 kgf), Lumbar: 49–588 N (5–60 kgf)
Treatment time	1–99 minutes (1 minute increments)
Safety function	Patient stop switch, Over-traction protection feature

Made in Japan 

Accessories

Accessories	Product name/Product number/Minimum ordering unit	
Blanket	OL5000-001	1
Cervical Traction Hanger	OL6600-004	1
Cervical Harness	OL6600-012	1
Pillow Heater	OL6600-003	1
Pillow (Gray)	OL6600-002	1
Side Pad (Slim) (Gray)	OL6600-011	1

- The colors in this brochure may not necessarily be the same as in the actual model due to printing conditions.
- Specifications and appearance are subject to change for improvement without notification.
- The content of this catalog is current as of April, 2023.

OG WELLNESS TECHNOLOGIES CO., LTD.

<http://www.og-wellness.com>

Head Office 1835-7, Miyoshi, Naka-ku, Okayama-shi, OKAYAMA 703-8261 JAPAN

Tokyo Head Office Kasumigaseki Building 4F, 3-2-5, Kasumigaseki, Chiyoda-ku, TOKYO 100-6004 JAPAN

Oku Factory 77, Mukoyama, Oku-cho, Setouchi-shi, OKAYAMA 701-4233 JAPAN